


This guide explains how to convert LabChart files into .txt files automatically. Before beginning, make sure that all LabChart files have been saved to a folder in a location you know (it's ok if the files are mixed in with TIFF files or something else in the folder, only the LabChart files will be converted).


### STEP 1:

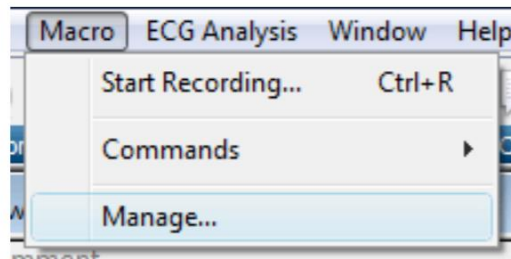

*Manage Macro from the menu bar.*

- In LabChart, select **Macro > Manage...** to bring up the Manage Macro dialog.


### STEP 2:

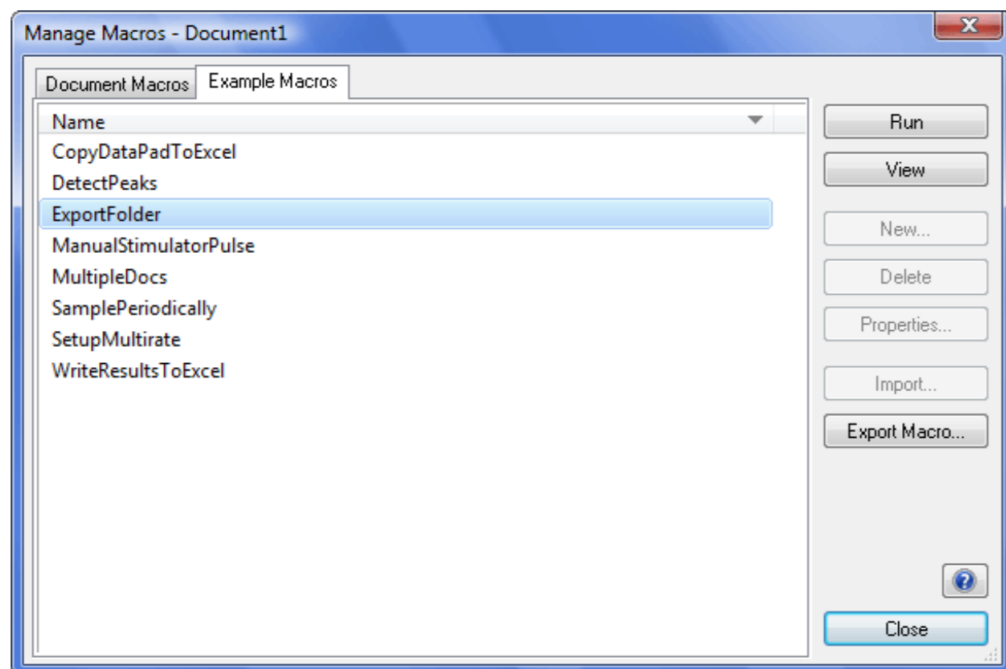

*ExportFolder Macro from the Manage Macros dialog.*

- Click on **Example Macros** tab.
- Select **ExportFolder** macro and press **View** to bring up the **ExportFolder** macro dialog.


### STEP 3:

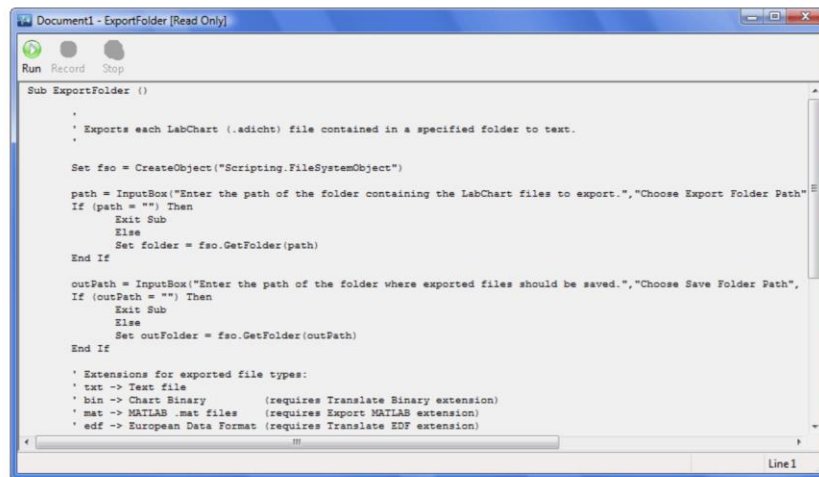

*ExportFolder Macro Dialog with Scripts.*


- Click **Run**.

### STEP 4:

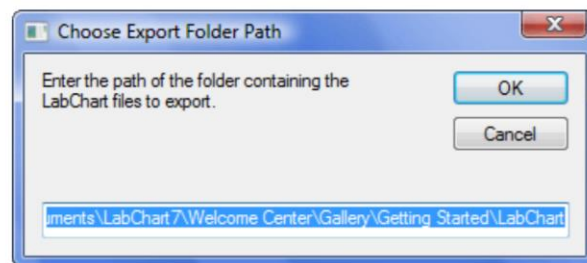

*Choose ExportFolder Path Dialog.*

- In the **Choose Export Folder Path dialog**, specify the folder containing files to be exported. Click **OK**.

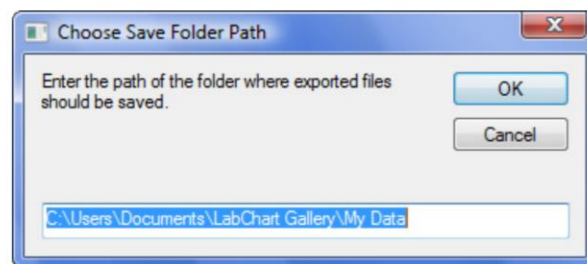

*Choose Save Folder Path Dialog*

- In the **Choose Save Folder Path dialog**, specify the folder where exported files are saved (the default location is the **My Data** folder in user's Documents folder). Click **OK**.

The **EXPORT FOLDER** is the folder where the original LabChart files are saved, and the **SAVE FOLDER** is any folder you want to save the .txt files into. It is this **SAVE FOLDER** that we will use for the next data analysis.